

AD-A178 877

DETERMINATION OF PHOSPHORUS IN WASTE WATER BY
INDUCTIVELY COUPLED PLASMA-ATOMIC EMISSION SPECTROMETRY
(U) NAVAL AIR DEVELOPMENT CENTER WARMINSTER PA A VAWMA
16 JAN 86 NADC-86138-68

1/1


UNCLASSIFIED

F/G 7/4

NL



(ND)
1 3/4
DTC

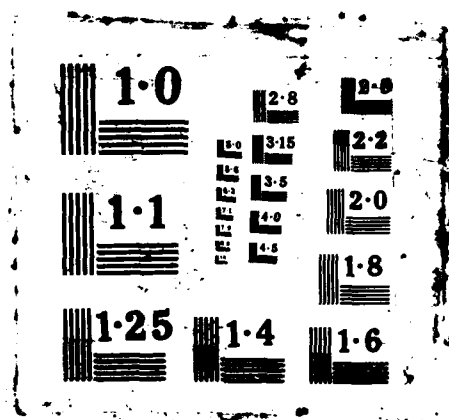

AD-A178 077

12

REPORT NO. NADC-86138-60



DETERMINATION OF PHOSPHORUS IN WASTE WATER BY INDUCTIVELY COUPLED PLASMA-ATOMIC EMISSION SPECTROMETRY

Asha Varma, PhD
Aircraft and Crew Systems Technology Directorate
NAVAL AIR DEVELOPMENT CENTER
Warminster, Pennsylvania 18974-5000

16 JANUARY 1986

FINAL REPORT
AIRTASK NO. 6201000001
Project No. 201000001A
Work Unit No. 233141
Project Element NO. APN

DTIC
ELECTE
MAR 19 1987
S D

Approved for Public Release; Distribution Unlimited

Prepared for
Office of Naval Technology
North Quincy Street
Arlington, Virginia 22217

AD-A178 077

87 3 18 037

UNCLASSIFIED

SECURITY CLASSIFICATION OF THIS PAGE

REPORT DOCUMENTATION PAGE

1a. REPORT SECURITY CLASSIFICATION UNCLASSIFIED			1b. RESTRICTIVE MARKINGS		
2a. SECURITY CLASSIFICATION AUTHORITY N/A			3. DISTRIBUTION/AVAILABILITY OF REPORT Approved for Public Release; Distribution is Limited.		
2b. DECLASSIFICATION/DOWNGRADING SCHEDULE N/A					
4. PERFORMING ORGANIZATION REPORT NUMBER(S) NADC-86138-60			5. MONITORING ORGANIZATION REPORT NUMBER(S) N/A		
6a. NAME OF PERFORMING ORGANIZATION Naval Air Development Center		6b. OFFICE SYMBOL (If applicable) 6062	7a. NAME OF MONITORING ORGANIZATION Office of Naval Technology		
6c. ADDRESS (City, State, and ZIP Code) Warminster, PA. 18974			7b. ADDRESS (City, State, and ZIP Code) North Quincy Street Arlington, VA 22217		
8a. NAME OF FUNDING / SPONSORING ORGANIZATION Office of Naval Technology		8b. OFFICE SYMBOL (If applicable) 225	9. PROCUREMENT INSTRUMENT IDENTIFICATION NUMBER		
8c. ADDRESS (City, State, and ZIP Code) North Quincy Street Arlington, VA 22217			10. SOURCE OF FUNDING NUMBERS		
			PROGRAM ELEMENT NO. APN	PROJECT NO. 201000001A	TASK NO. 6201000001
			WORK UNIT ACCESSION NO. 233141		
11. TITLE (Include Security Classification) DETERMINATION OF PHOSPHORUS IN WASTE WATER BY INDUCTIVELY COUPLED PLASMA ATOMIC EMISSION SPECTROSCOPY					
12. PERSONAL AUTHOR(S) Asha Varma					
13a. TYPE OF REPORT Final		13b. TIME COVERED FROM Nov 85 to Dec. 85		14. DATE OF REPORT (Year, Month, Day) Jan. 16, 1986	
15. PAGE COUNT					
16. SUPPLEMENTARY NOTATION N/A					
17. COSATI CODES			18. SUBJECT TERMS (Continue on reverse if necessary and identify by block number)		
FIELD	GROUP	SUB-GROUP	Inductively coupled plasma-atomic emission spectroscopy (ICP-AES), Waste, Standard addition		
19. ABSTRACT (Continue on reverse if necessary and identify by block number) Under the EPA regulations, waste water should not contain more than 5-6mg/L of total phosphorus before it is discharged into a public sewer system. The conventional methods for determining phosphorus in waste are time consuming and also require pre-treatment of samples to convert all the phosphorus to the ortho-phosphate form for analysis. Inductively coupled plasma-atomic emission spectrometric method (ICP-AES) has been successfully used to determine phosphorus in waste. This method is fast, accurate and does not require sample pretreatment prior to analysis. <i>Keywords</i>					
20. DISTRIBUTION/AVAILABILITY OF ABSTRACT <input checked="" type="checkbox"/> UNCLASSIFIED/UNLIMITED <input type="checkbox"/> SAME AS RPT. <input type="checkbox"/> DTIC USERS			21. ABSTRACT SECURITY CLASSIFICATION UNCLASSIFIED		
22a. NAME OF RESPONSIBLE INDIVIDUAL Dr. Asha Varma			22b. TELEPHONE (Include Area Code) (215) 441-3975		22c. OFFICE SYMBOL CODE 6062

DD FORM 1473, 84 MAR

83 APR edition may be used until exhausted.

All other editions are obsolete

SECURITY CLASSIFICATION OF THIS PAGE

U.S. Government Printing Office: 1988-639-012

0102-LF-014-6602

UNCLASSIFIED

TABLE OF CONTENTS


Section	Page
INTRODUCTION	1
EXPERIMENTAL PROCEDURE	2
RESULTS	2
CONCLUSION	3
REFERENCES	10

LIST OF TABLES

Table No.		Page
1.	Instrumental set-up parameters for phosphorus analysis	4
2.	Phosphorus in waste water	5
3.	Phosphorus in waste water by standard addition method	6
4.	Comparative values for phosphorus	7

LIST OF FIGURES

Figure No.		Page
1.	Graphic spectrum of phosphorus 200, 100 and 20 mg/L	8
2.	Graphic spectrum of phosphorus 20, 10 and 5 mg/L	9

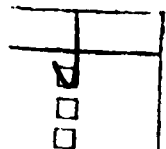

Di t

A-1

Codes

INTRODUCTION

Phosphorus is one of the important elements to be detected in environmental, biological, geochemical and metallurgical products. Murphy and Riley¹ proposed an application of phospho-molybdenum heteropoly method for the determination of phosphorus in sea water, and it became one of the most popular analytical techniques. The standard methods for phosphorus²⁻⁴ determination are time consuming and require pretreatment of samples for conversion of all the phosphates to the ortho-phosphate form for analysis. Although either direct or indirect atomic absorption spectrometric method⁵⁻¹⁰ is a fast method of analysis, the phosphorus concentration in some samples is too low to be detected by this method.

Inductively coupled plasma-atomic emission spectrometry (ICP-AES) method is characterized by a high detection efficiency and good reproducibility, and is an accepted technique for quantitative analysis. This method provides speed, accuracy and reliability for the determination of phosphorus in steel¹¹, rocks¹², water^{5,13,14}, waste¹⁵ etc. Wings et al¹³ reported two wavelengths (213.618nm and 214.914nm) for phosphorus determination using ICP-AES method. It has been reported¹⁶ that a partially resolved copper emission line at 214.90nm interferes with the 214.914nm phosphorus line. Other methods using ICP detector¹⁷ and hyphenated ICP techniques¹⁸⁻²¹ for phosphorus determination are also reported in the literature.

In this report ICP-AES method for the determination of phosphorus in waste water is presented. Phosphorus in 1-5 mg/L quantity can be determined with $\pm 1.2\%$ error.

EXPERIMENTAL PROCEDURE

Reagents:

Phosphorus standard	1000 mg/L
Phenolphthalein indicator	
Sodium hydroxide	5.0 N
Sulfuric acid	10.0 N

Samples:

Sample solutions were identified as,

Raw sewage(RS)
Raw industrial waste (RIW)
Final effluent (FE)
Treated industrial waste (TIW)

Equipment

The instrument used for all the measurements was a Perkin-Elmer ICP/6000 spectrometer equipped with a 7500 computer. The instrument set-up parameters for the phosphorus analysis were as described in Table 1. The analysis was performed with and without purging the ICP system.

Procedure:

Standard phosphorus solution was prepared from a stock solution of phosphorus (Spex Industries) by serial dilution with deionized distilled water. Deionized distilled water was also used as a blank and a rinse after each sample analysis to minimize the memory effect. Use of glassware was avoided as much as possible, if and when used it was rinsed with diluted hydrochloric acid.

Since the amount of phosphorus presented in the sample was very low as determined from the preliminary runs, a 20 mg/L phosphorus standard solution was used for all the analyses. Sample solutions were analysed in the as-received condition. Any solid residue if present in the sample solutions was removed by filtration or it would clog the nebulizer. It was found that the addition of a few ml of acid helped in dissolving the suspended solid. Colloidal solutions were easy to analyse without any interference. 10-15 sets of each sample were analyzed for consistent values. A standard addition method was used to make sure that the amount of phosphorus determined was correct. For comparison purposes samples were also subjected to pretreated as suggested under the Standard Methods for Phosphorus and then analyzed under the similar conditions.

RESULTS

The concentration of phosphorus found in raw and treated sewage water is reported in Table 2. An accurate determination of phosphorus was also made by the standard addition method and data obtained are presented in Table 3. The comparative values for phosphorus present in the sewage water by the direct analysis and the standard addition methods are reported in Table 4. The phosphorus concentration found in the treated and untreated samples was exactly the same.

A typical graphic spectrum is shown in Fig. 1 for 200, 100 and 20 mg/L and in Fig. 2 for 20, 10 and 5 mg/L phosphorus concentrations. A peak left of the standard phosphorus peak is identified as also due to phosphorus (identified as #5 in Figs. 1 & 2). Another peak (identified as #6 in Figs. 1 & 2) was observed to the right of the phosphorus peak and was identified by the Computer Search as due to As, Cu, Fe, Ni, V and Zn. No interferences were observed from any of these elements in the determination of phosphorus by this method.

CONCLUSION

Results obtained from the ICP-AES analysis of sewage water by direct analysis and the standard addition methods prove that this method can be used for routine, fast and accurate analysis of phosphorus. A quantitative value of 1 mg/L phosphorus can be obtained with $\pm 1.2\%$ error.

Table 1-Instrumental Set-Up Parameters For Phosphorus Analysis

Plasma	Argon, high purity
Plasma height	15 nm
Plasma flow	12 L/min.
Auxillary flow	1 L/min.
Nebulizer	24-25 psi
RF generator	27.12 MHz
Incident power	1250
Reflectance	5 w
Purge	Nitrogen or argon gas
Wavelength	213.618 nm
Range	1.0 nm
Gain	Automatic
Background corr.	Both sides of the peak
Detection limits	0.076 mg/L at 213.618 nm
	0.05 mg/L at 214.914 nm

Table 2 - Phosphorus in waste water

Sample name	Concentration of phosphorus found, mg/L	Standard deviation
Raw sewage (RS)	8.34 8.15 8.91	0.28 0.05 0.51
Final effluent (FE)	5.56 5.42 5.56	0.34 0.24 0.23
Raw industrial waste (RIW)	0.52 0.50 0.52 0.53	0.09 0.14 0.13 0.21
Treated industrial waste (TIW)	0.00	0.00

Table 3 - Phosphorus In Waste Water By Standard Addition Method

Sample name	Standard P added, mg/L	Total P found, mg/L	Actual P present mg/L	Standard deviation
RS	1.0	9.809	8.809	0.2
	1.0	9.589	8.589	0.5
	2.0	10.606	8.606	0.5
	2.0	10.525	8.525	0.5
	3.0	11.86	8.53	0.3
FE	2.0	7.53	5.53	0.2
	2.0	7.59	5.59	0.2
	2.25	7.84	5.59	0.4
	2.50	8.05	5.55	0.4
RIW	2.5	3.12	0.62	0.3
	4.0	4.792	0.79	0.2
	4.0	4.786	0.79	0.4
	5.0	5.634	0.63	0.4
	5.0	5.746	0.63	0.3
TIW	2.0	2.0	0.00	0.0
	2.0	1.8	-0.02	0.2
	2.5	2.5	0.00	0.0
	2.5	2.52	+0.02	0.2
	3.33	3.26	0.07	0.1

Table 4 - Comparative Values For Phosphorus

Sample name	Phosphorus in mg/L		Difference \pm
	Direct analysis	Standard addition	
RS	8.67	8.61	-0.06
FE	5.53	5.56	+0.03
RIW	0.64	0.69	+0.05
TIW	0.070	0.00	± 0.00

Graphics Spectrum Mode

85/06/24

09:30

Spectra Name: PHOSPHORUS Wavelength: 213.618 nm

Range: 1.00 nm

Gain: 538

Background Intervals: -0.045 nm + 0.035 nm

Read Delay: 20

Remarks: Standard 200, 100 and 20 mg/L phosphorus solution

1 Standard 200

2 Blank

3 100 mg/L

4 20 mg/L

5

6

7

8

9

10

11

12

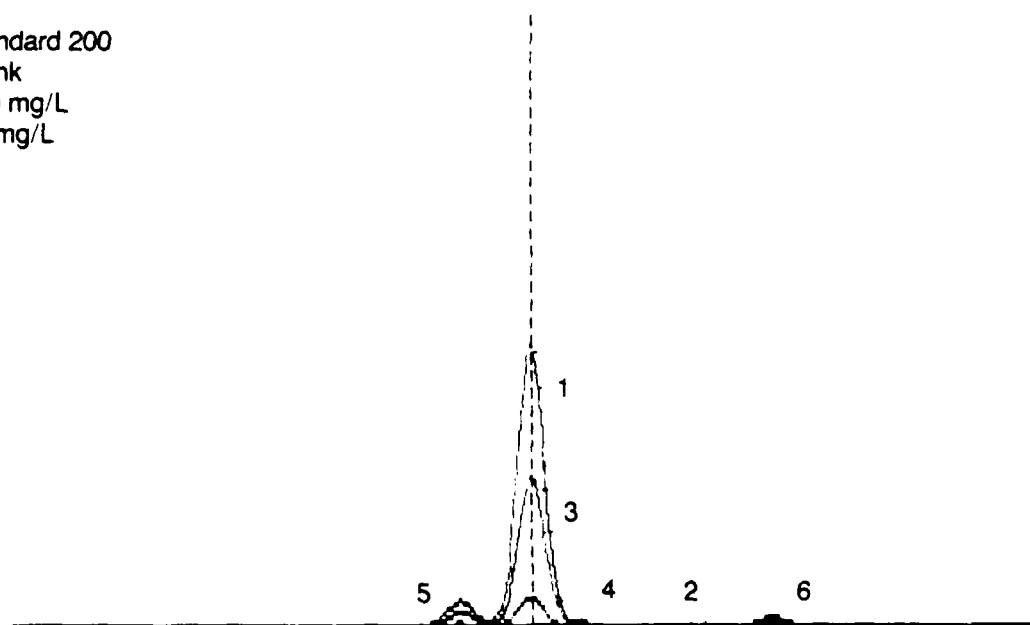


Fig. 1. Graphic Spectrum of Phosphorus, 200, 100 and 20 mg/L

Graphics Spectrum Mode

85/06/24

09:46

Spectra Name: PHOSPHORUS Wavelength: 213.618 nm Range: 1:00 nm Gain: 638
Background Intervals: -0.045 nm + 0.066 nm Read Delay: 20
Remarks: Standard phosphorus in 20, 10 and 5 mg/L concentration

- 1 Standard 20ppm
- 2 Blank
- 3 10ppm
- 4 5ppm
- 5
- 6
- 7
- 8
- 9
- 10
- 11
- 12

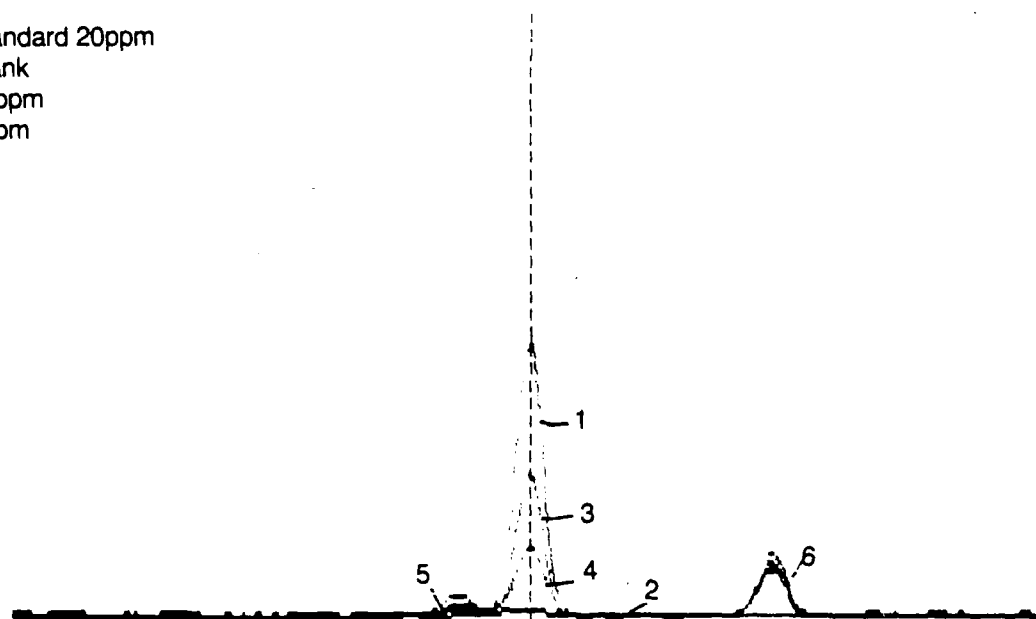


Fig. 2. Graphic Spectrum of Phosphorus, 20, 10 and 5 mg/L

REFERENCES

1. Murphy, J. and Riley, J., Anal. Chim. Acta 27, 481, 624(1975).
2. APHA Standard Methods, 14th Edition, pp473, 481, 624 (1975).
3. Annual Book of ASTM Standards, Part 31, "Water", Standard D515-72, p388 (1976).
4. Environmental Protection Agency, Guidelines establishing Test procedures for the Analysis of Pollutants; Proposed Regulations, Fed. Regist. 03 Dec. 1979, 44(233), 68424-69573; ICP Inf. Newsl. 5, 528 (1980).
5. Pongar, R. W., and Thomson, G.R., Jarrell Ash Plasma Newsletter 1 (2), 5, April 1978.
6. Vigler, M. S., Strecker, A. and Varnes, A., Appl. Spectrosc. 32 60-62 (1978).
7. Driscoll, D. J., Clay, D. A., Rogers, C. H., Jungers, R. H. and Butler, F. E., Anal. Chem. 50, 767-79 (1978).
8. Kubota, T., Ueda, T. and Okutani, T., Bunseki Kagaku 33, 633 (1984).
9. Grease, F. and Mardi, J. G. Analusis 12, 467 (1984).
10. Lin, S. W. and Julshamn, K., Anal. Chim. Acta 158, 199 (1984).
11. Wallace, G. F., Hoult, D. W. and Ediger, R. D., At. Spectrosc. 1(4), 120 (1980).
12. Hirata, S., Bunseki Kagaku 33, T64 (1984).
13. Winge, R. K., Peterson, V. J. and Fassel, V. A., Appl. Spectrosc. 33(3), 206 (1979).
14. Floyd, M. A., Haloums, A. A., Morrow R. W. and Farrar, R. B., Am. Lab. p84, March 1985.
15. Heinrich, D. L., Lab Equipment p42, OCT. 1983.
16. Nygaard, D. D., Chase, D. S., Leighty, D. A. and Smith, S. B., Anal. Chem 56, 424 (1984).
17. Heine, D. R., Denton, M. B. and Schlabach, T D., Anal. Chem. 54, 81 (1981).
18. Matsumoto, K. and Fuwa, K., Bunseki Kagaku 30, 188 (1981).
19. Yoshida, K., Haraguchi, H. and Fuwa, K., Anal. Chem 55, 1009 (1983).
20. Urasa, I. T., Anal. Chem. 56, 904 (1984).
21. Hashimoto, S., Fujiwara, K. and Fuwa, K., Anal. Chem. 57, 1305 (1985).

END

4-87

DTIC